



Mecca Care Shower Seat 960×330mm

Nero

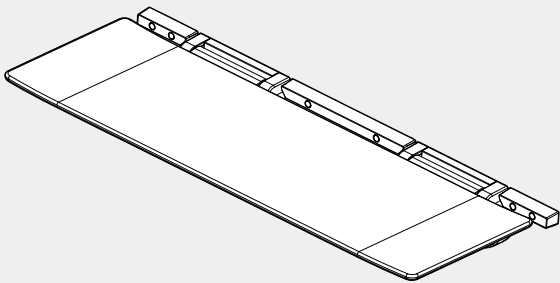

Version: 202606

Drawing

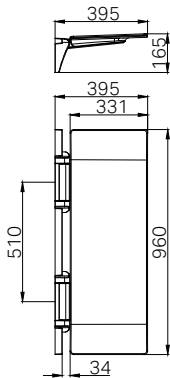

Features

1. Ergonomic design folds neatly up against wall.
2. High quality, high strength materials and finishes.
3. Compact size suitable for all shower sizes.
4. 150kg capacity when installed into brick/masonry or multiple wall studs.
5. AS1428.1 Compliant for DDA bathrooms.
6. Gun metal, brushed gold and brushed bronze colours are PVD coated for long term durability.

Finish & SKU

Chrome
NRCR0002CH

Gun Metal
NRCR0002GM

Brushed Bronze
NRCR0002BZ

Matte Black
NRCR0002MB

Brushed Gold
NRCR0002BG

Brushed Nickel
NRCR0002BN

Warning


These rails are designed to be mounted onto wall studs or masonry walls. Installation is recommended to be completed by qualified trades person who complies with the AS 1428.1-2009 design for access and mobility standard.
It is recommended rail location is determined by a qualified occupational therapist or trained professional to the AS 1428.1-2009 standard. Correct installation is required to ensure safe, long term operation. Incorrect mounting can result in reduced capacity of the rail which can result in injury to user and damage to the wall.

*Dimensions are nominal measurements only.

Installation Instruction


Cautions

- 1: Please confirm whether the wall is a load-bearing before drilling;
- 2: Before drilling, please confirm whether there is a water pipe at the hole position;
- 3: Check whether the product accessories are complete before installation;
- 4: Make sure the drill bit is 10mm and 6mm in size.
- 5: Make sure the hole depth is 6.5cm and 3.5cm.
- 6: The recommended height from floor to top of seat is about 470-480mm.

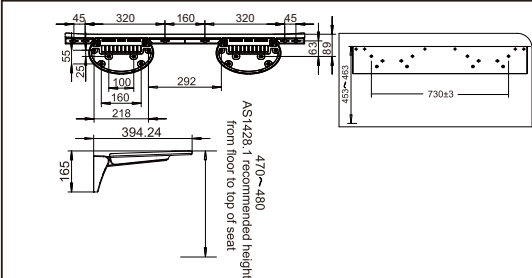


1. Ensure wall is suitable for installing the shower seat before you start installation.
Place wall fixing plate on wall in correct position and mark out holes.


Tips: Recommend height from floor to top of seat is 470-480mm. Then recommend height from floor to wall fixing plate is 453-463mm, refer to left picture.

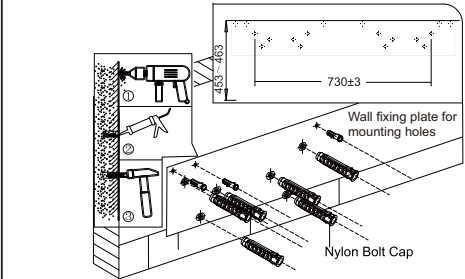


2-1. Drill 12 holes on the wall based on the screw holes of our shower seat (with a drill bit of 10mm and a hole depth of 6.5cm).
2-2. Drill 6 holes on the wall based on the screw holes of our shower seat (with a drill bit of 6mm and a hole depth of 3.5cm). Apply silicone to the hole and use a hammer to knock the Nylon bolt cap into the hole on the wall.

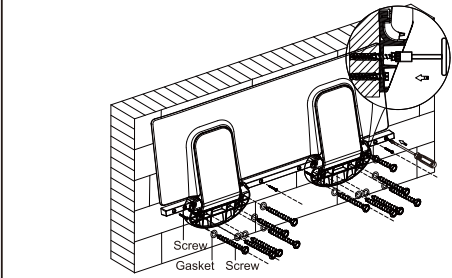

3. Place the shower seat on the wall, align the screw holes with the wall holes, and then tighten the screws with a socket wrench.

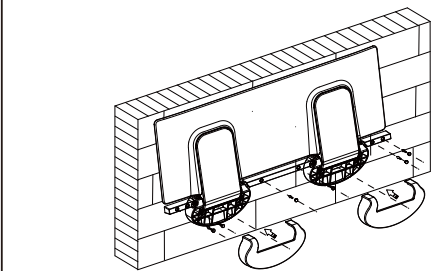

4. Tighten the base cover.

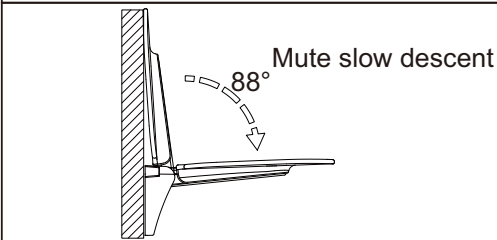

*Note: you can choose the style without slow down function, and the installation method is the same.

5. After installation, the seat shall fold in an upwards direction and when folded the grabrail shall be accessible. Unfold it When in use.

*Note: 2-1. The drilling bit is 10mm, and the drilling depth is greater than or equal to 6.5cm
2-2. The drilling bit is 6mm, and the drilling depth is greater than or equal to 3.5cm

Nero

